CI/SfB			
	(74.26)	Х	(U35)




www.chilterninvadex.co.uk

Showering and Bathing

Level, Easy and Step-in **Shower Range**


Easy bathing solutions for the home, care environment and hospital

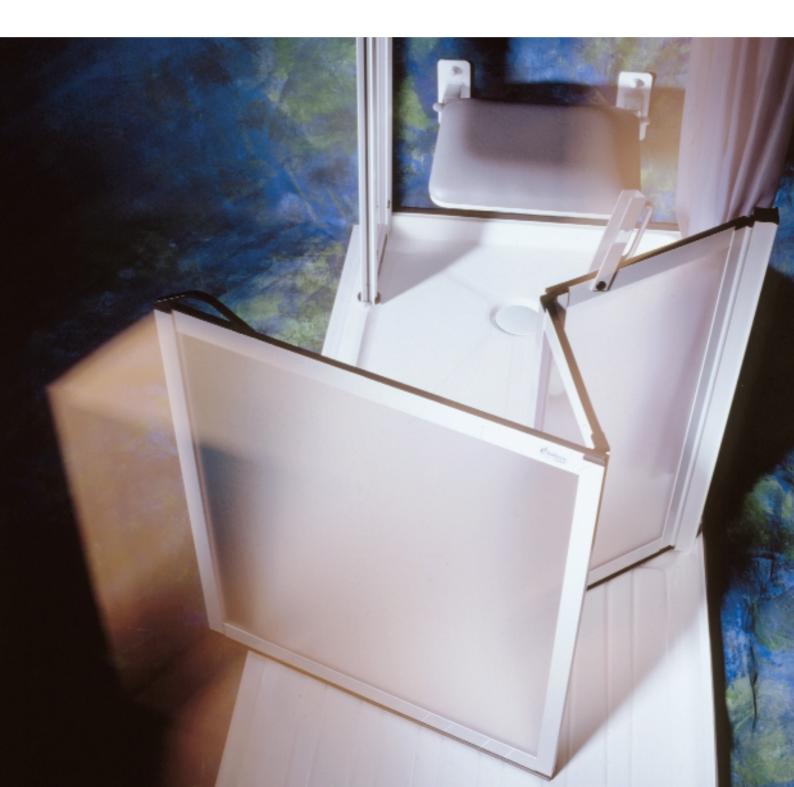

Shower Tray Selector

To aid selection of the most appropriate shower tray to suit the floor construction, user needs and space considerations please use this easy to follow chart. Features are denoted by coloured shape.

1. Drainage Options

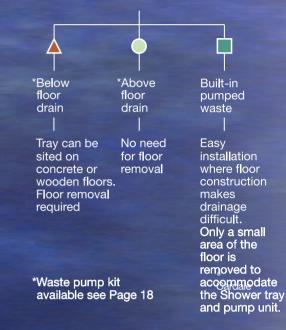

2. What are clients/family needs (now and in the future)?



- * Subject to type (attendant or self propelled model) and width of shower chair chosen. See page 23 for shower chair compatibility.
- 3. What will fit the site taking into consideration answers from above? (detailed dimensions are contained on page 20-22)

Select shower tray

Shower Tray

Sandwell MK3 655/760 Westdale 655/760/1000 Abbeydale 760/850/950/1050 Dovedale MK2 Cardale MK2 Edale MK2 Berriedale Ferndale



athing Showering and

On the reverse of this flap is a Shower Tray Selector chart, designed to be left unfolded as a quick guide for easy reference.

Level, Easy and Step-in Shower Range

Solutions for Living

Chiltern Invadex provide a wide range of products designed to maintain independence, aid mobility and improve the quality of daily life.

Product and service quality

All shower trays and enclosures are manufactured from high quality durable materials, which retain their appearance, even with regular use. Shower trays are guaranteed against breakage for 5 years, while virtually unbreakable Invaglaze panels are guaranteed against yellowing for 5 years.

Customer service centres can provide technical support, query resolution and arrange showroom bookings.

Survey and Installation

A full installation service is available for all Chiltern Invadex Shower Trays and Enclosures. A detailed site survey and scale drawing can be provided without charge by a Technical Sales Representative. Please call your nearest Chiltern Invadex Office for details.



Showering and Bathing

Level, Easy and Step-in **Shower Range**

Solutions for living



LEVEL ACCESS P 2-1

The Level access range of showers combine stylish and innovative design with floor level access, without the need to overcome ramps or steps. Simple installation on most concrete or wooden floors.

Sandwell

Ideal bath replacement with floor level access for ambulant or shower chair users. Incorporates a showering and drying area

Westdale

Level access for ambulant or shower chair users with size options to suit different size bathrooms

Cardale MK2

Easily installed on virtually any surface and for ambulant or shower chair users

ACCESS P 8-13 Е ASY

The Easy access range of showers combine stylish looks with ease of installation and includes models for the smaller bathroom, for wheeled shower chair and for walk in access.

- Abbeydale
- Dovedale MK2 Ideal for ambulant users
- for ambulant or shower chair users Edale MK2



Choice of sizes with extra space

<u>STEP-IN SHOWERS P14-16</u>

Designed for those able to step into the shower, this versatile range of shower trays is available with a choice of drainage positions and trav dimensions, creating a flexible solution to suit variations in bathroom layout.

Berriedale

Ideal for step-in users, simple to install and with extra length

Ferndale

Ideal for the bathroom with limited space, and with a shallow lip for step-in users



ACCESSORIES P 17-18

Units and Grab rails.

TECHNICAL INFORMATION P 19-23

Technical information on all products is clearly set out for easy specification on pages 19-23.









Sandwell MK3 655/760



S

An ideal bath replacement with floor level access for ambulant or shower chair users

Combining floor level access with stylish design the Sandwell MK3 has no ramps or steps for the bather to overcome and incorporates a showering and drying area.



Tri-Fold Doors

Please take advantage of our free on site surveying service.

Package includes

- Shower tray (Left or right hand)Door enclosure (see page 20
- for options)
 Tray edging strip
- Inva-waste water trap
- Curtain and curtain rail
- Assistance pole

Features

- *Level Access* With no ramp or steps to overcome and ideal for walk-in or wheeled shower chair users.
- Choice of width See page 20 for options.
- *Simple installation -* Fits most wooden or concrete floors.



С

Slide and Fold Doors

- *Tray colour matching choice* A standard range available at no extra cost.
- Easy clean waste trap.
- *Trimming the tray* Length can be easily cut to fit the space.
- Choice of left or right hand drain options.
- *Easy grip assistance pole -*Assistance pole to curtain rail is provided as standard to help get in and out of the shower (option to the ceiling also available).
- Textured floor surface for added safety.
- *Tray edging strip* Provides a neat watertight join between tray and flooring.
- Door options -
 - Tri-fold door arrangement.
 Slide and fold doors, which cover the length of the tray and drying area (above right).







Self Propelled


Shower Chair



Pumped waste option see page 18

Static Shower Seat

Attendant Propelled Shower Chair

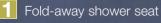

Grab rails 2

3

Choice of shower units from a range of leading manufacturers

Sandwell MK3 655/760

A

An ideal bath replacement with floor level access for ambulant or shower chair users

3

1

Westdale 655, 760, 1000

Level access for ambulant or shower chair users

Different sized bathrooms with limited space for mobility are easily overcome with the Westdale range of level access showers (with no ramps or steps to negotiate) available in three alternative sizes.



Please take advantage of our free on site surveying service.

Package includes

- Shower tray (Left or right hand)
- Tri-fold door enclosure
- Tray edging strip
- Inva-waste water trap
- Curtain and curtain rail
- Assistance pole

Features

- Level Access Floor level access with no ramp or steps to overcome and Ideal for walk-in or wheeled shower chair users.
- Simple installation Fits most wooden or concrete floors.



С

S

- Choice of three sizes Options to suit all bathroom layouts. See page 20 for options.
- *Tray colour matching choice* A standard range available at no extra cost.
- Easy clean waste trap.
- A choice of drain options -Available in left or right hand drain options.
- *Easy grip assistance pole -*Assistance pole to curtain rail is provided as standard to help get in and out of the shower (option to ceiling also available).
- Textured floor surface provides added safety.
- *Tray edging strip* Provides a neat join between tray and flooring.







Self Propelled

Shower Chair



Pumped waste option see page 18

Static Shower Seat

Attendant Propelled Shower Chair

Fold-away shower seat

2 Grab rail

3

Choice of shower units from a range of leading manufacturers

[]

Westdale 655, 760, 1000 Level access for ambulant or shower chair users

Res

3

1.1

Cardale MK2 with pumped waste



S

Easily installed almost anywhere with level access for ambulant and shower chair users

Installing a shower tray on a solid floor presents no problem for the Cardale MK2, which has been innovatively designed for installation on almost any surface, or in any room, with the minimum of structural disruption. The Cardale MK2 also offers level access, with no steps or ramps to overcome.



Full Height Panel Option

Please take advantage of our free on site surveying service.

Package includes

- Shower tray (Left or right hand)
- Tri-fold door enclosure
- Pump, pump cover & wiring block
- Curtain and curtain rail
- Assistance pole

Features

• Simple installation - fits most wooden or concrete floors.



С

Half Height Panel Option

- Integral pumped waste system -The unique pumped waste system (with quiet 12 volt pump) overcomes difficult drainage problems.
- Choice of installation method -The Cardale MK2 can be supplied without panels, with half, or full height wall panels. Panel options, reduce installation time, and conceal pipe work.
- *Tray colour matching choice* A standard range available at no extra cost.
- Half height door enclosures -Carer friendly door height allows a carer to assist without getting wet.
- *Easy grip assistance pole -*Assistance pole to curtain rail is provided as standard to help get in and out of the shower (option to ceiling also available).
- Textured floor surface for added safety.







Static Shower Seat

Attendant Propelled Shower Chair Self Propelled Shower Chair

1 Fold-away shower seat

2 Grab rails

Choice of shower units from a range of leading manufacturers

0

Cardale Level access for ambulant and shower chair users

3

64

1

6

Abbeydale 760, 850, 950 & 1050 Choice of sizes with easy access

С

С

F

S

Available in four sizes the Abbeydale range is a shower tray and enclosure with extra room to accommodate both wheeled shower chair and walk-in users. Low entry access makes getting in and out of the shower simple.

S



Please take advantage of our free on site surveying service.

Package includes

- Shower tray (Left or right hand - 760 model only)
- Tri-fold door enclosure (Left or right handed - see page 21 for options)
- Inva-waste water trap
- Tray edging strip
- Curtain and curtain rail
- Assistance pole

Features

- Choice of four sizes -See page 21 for options.
- Simple installation The Abbeydale is quick and easy to install.







Static Shower Seat

Attendant Propelled Shower Chair Self Propelled Shower Chair

Pumped waste option see page 18

- *Easy access* Easy low entry access and neat tray edging strip provide an unobtrusive join between floor and shower tray.
- Below floor waste Tray sits on the floor (only a small area of the floor is removed for the drain).
- Door options Tri-fold doors are standard, hinged on the left or right.
- Half height door enclosures -Carer friendly door height allows a carer to assist without getting wet.
- Textured floor surface for added safety.
- Easy grip assistance pole -Assistance pole to curtain rail is provided as standard on the 760 model (two full height assistance poles on the 850, 950 and 1050 models) to help the user get in and out of the shower (option to ceiling also available).
- Tray colour matching choice -A standard range available at no extra cost.

8





Fold-away shower seat

2 Grab rails

Choice of shower units from a range of leading manufacturers

Abbeydale 760, 850, 950 & 1050 Choice of sizes with easy access

Dovedale MK2

S



S

Easy access for the ambulant user

Compact in size, yet simple to install the Dovedale MK2 features a removable easy access ramp and a host of user options to help make the shower more accessible.



Please take advantage of our free on site surveying service.

Package includes

- Shower tray
- Door enclosure (see page 21 for options)
- Removable ramp
 Curtain and curtain rail*
- *not supplied with full height panels • Assistance pole

Features

- *Simple installation -* The Dovedale MK2 is quick and easy to install.
- Easy access Incorporates a removable easy access ramp.



С

С

- Above floor waste No need to remove the floor, as the waste is discharged above floor with 4 optional drainage positions.
- *Door options* Single or Tri-fold doors are optional.
- *Half height door enclosures* -Carer friendly door height allows a carer to assist without getting wet.
- *Panel options* Half height side panels supplied as standard. Full height side panels to order.
- Textured floor surface for added safety.
- *Easy grip assistance pole -*Assistance pole to curtain rail is provided as standard to help get in and out of the shower (option to ceiling also available).



Static Shower Seat

Pumped waste option see page 18



Fold-away shower seat

2 Grab rails

Choice of shower units from a range of leading manufacturers

Dovedale MK2

Easy access for the ambulant user

Edale MK2

S



The Edale MK2 sits neatly on the floor surface without the need to remove a section of the floor for drainage, making installation simple. The low entry detachable ramp (optional) is carefully graduated to help make the shower accessible for wheeled shower chair or walk-in users.



Please take advantage of our free on site surveying service.

Package includes

- Shower tray (Left or right hand)
- Tri-fold door enclosure
- Assistance poles
- Curtain and curtain rail
- Removable ramp is optional

Features

• *Easy to install with above floor waste* - No need to remove the floor, as the waste water is discharged above the floor, which makes plumbing simple.



С

S

С

 Adjustable and safe - The horizontal hand rail can be positioned at a comfortable height to suit the user, providing assistance when entering or sitting.

Assistance pole to curtain rail is provided as standard to help get in and out of the shower (option to the ceiling also available).

- *Door options -* Tri-fold doors supplied as standard.
- Half height door enclosures -Carer friendly door height allows a carer to assist without getting wet.
- *Tray colour matching choice* A standard range available at no extra cost.
- Textured floor surface for added safety.







Self Propelled

Shower Chair



Pumped waste option see page 18

Static Shower Seat

Attendant Propelled Shower Chair





Fold-away shower seat

2 Grab rails

3 Choice of shower units from a range of leading manufacturers

1

Ideal for any surface and for walk-in or wheeled shower chair users

Edale MK2

3

2

Berriedale

Ε

Т

Ρ

S



Е

R

Ideal for step-in users and with extra leg room

Ν

S

Η

0

W

Ideal for the step-in user with its low access, the Berriedale is simple to install, offering a practical solution for people who need extra leg room. The tray sits neatly on the floor surface making installation quick and easy.



Please take advantage of our free on site surveying service.

Package includes

- Shower tray (Left or Right hand)
- Door enclosure
- Assistance pole
- Curtain and curtain rail

Features

- Low access Ideal for step in users.
- *Side drain outlet -* Easily accommodates variations in bathroom layout.



- *Simple installation* The Berriedale sits neatly on the floor surface.
- Half height door enclosures to order Ideal for assisted bathing.
- Handed drainage Left and right hand drain options.
- *Tray colour matching choice* A standard range available at no extra cost.
- Textured floor surface for added safety.
- Choice of shower units -Available for use with this product.



Static Shower Seat

Pumped waste option see page 18

Grab rails

2 Choice of shower units from a range of leading manufacturers

2

Ke

1

10

Berriedale

Ideal for step-in users and with extra length

25



Ferndale



Ideal for step-in users and for the smaller bathroom

Ideal for the step-in user with its low access the Ferndale can replace a conventional bath and sits neatly on the floor surface, making installation quick and easy.



Package includes

Right hand)

door options)

Features

Assistance pole

Curtain and curtain rail

White shower tray (Left or

Door enclosure (Tri-fold or single

• Low access - access height of

• Side drain outlet - left or right option to accommodate variations

in bathroom layout.

120mm is ideal for step-in users.

• Simple installation - the Ferndale

sits neatly on the floor surface.



Single door option

- *Tray colour matching choice* a standard range available at no extra cost.
- Textured floor surface for added safety.
- Choice of shower units from a range of leading manufacturers.
- *Wall panel options* lightweight panel options reduce installation time and conceal pipework.
- *The Ferndale shower tray* can be supplied without an enclosure, with a curtain rail/curtain option.

NOTE: Illustration shows wall panel option. The Ferndale can be tiled without wall panels.





Static Shower Seat

Pumped waste option see page 18

Complementary Accessories

A full range of complementary accessories is available



TransAqua shower chairs

Combining durable stainless steel with modular construction the TransAqua range offers flexibility for both user and carer.

- Attendant and self propelling models
- 18" and 21" width options
- Height and width adjustable
- Soft padded seat completely sealed against water
- Suits standard WC's/Clos-o-Mat
 Easily dismantled without tools





Shower Seats

Available with flat-padded, horseshoe or aperture seat with anti-bacterial polypropylene base.

To see the complete range please ask for a copy of the Chiltern Invadex Shower Seats brochure.



Grab Rails

Our strong and secure grab handles and rails are designed to provide extra support in the bathroom.

To see the complete range please ask for a copy of the Chiltern Invadex Accessories brochure.



Shower and Mixer units

Choose your shower from a selection of leading manufacturers including Mira and Triton.

To see the complete range please ask for a copy of the Chiltern Invadex Accessories brochure.



New Wet floor kit

The new Wet Floor Kit features a quiet and filterless pump, which quietly pumps waste water away in a wet floor installation, where conventional drainage may be difficult.

Kit contains

- Flexible drain hose
- Wiring block
- Quiet 12 volt pump
- Stylish corner cover
- Clamping plate to grip floor covering.





Shower waste pump options available

Portascreen Shower Screen


Containing the spray from a shower is made simple with the Portascreen shower guard, which allows a carer to assist on an open floor without getting wet. Highly portable (8.5 kg) and easy to handle, the Portascreen folds flat for storage.



Warm air Body Dryer MK2

The Chiltern Invadex Body Dryer is ideal for those who find towel drying difficult or uncomfortable. Both simple to operate and hygienic to use, the Body Dryer soothes and relaxes aching joints and muscles.

- Slim profile makes the Body Dryer suitable for installation in any showering area or shower tray and enclosure.
- Tear drop pattern holes discharge warm, pressurised air that tones and invigorates the body.
- The touch control can be mounted in a convenient position, even on the floor, making it easy to use and helping to maintain independence for the user.

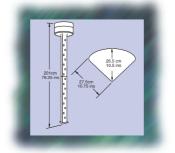

Shower waste pump

To overcome situations where access to drainage is difficult a shower waste pump allows waste water to be pumped away, without the need to rely on a gravity fall.

All Chiltern Invadex shower trays can be fitted with a Shower Waste Pump (except Cardale where it comes as standard).

- Quiet and filterless pump
- Simple to install and maintain
- Audible alarm to indicate blockages
- Wiring block included



Invaflo Shower Drain

Easy to install clamping plate for PVC floor membrane with built in gully for surface access. No earth bonding required under electrical regulations. Available in a choice of colours and lengths.



Technical specifications

Casing Electrical supply	
Power consumption Nominal output	
Average air outlet	

Safety cutouts installed

temperature

Fan

Moulded ABS 220/240 volt 50/60Hz 1PH 12.5 Amps 2.9 kilowatts

70° Centigrade 2 speed

TECHNICAL INFORMATION

Technical Information




How to work out whether you need a left or right handed shower tray

The position of the waste outlet determines the handing of the tray (as you look in on to the tray from the longest entry side) so if the waste outlet is on the right the tray is a right handed shower tray.

Space requirements


Contact your nearest Chiltern Invadex Office if you are unsure about space requirements for shower trays and enclosures.

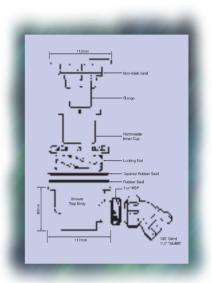

Inva Waste

	Important note <i>Power showers</i> Recommended maximum water delivery rate for Chiltern Invadex level, easy and step-in access shower enclosures.								
	Model	Flow Rate							
	Sandwell MK3	9lt/min (2 gal/min)							
	Westdale	9lt/min (2 gal/min)							
	Abbeydale	9lt/min (2 gal/min)							
	Cardale MK2	9lt/min (2 gal/min)							
	Dovedale MK2	9lt/min (2 gal/min)							
	Edale MK2	9lt/min (2 gal/min)							
	Berriedale	9lt/min (2 gal/min)							
	Ferndale	9lt/min (2 gal/min)							
L	EVEL	A C C E S S							
	Sandwell MK3 655, 760Cardale MK2	 Westdale 655, 760, 1000 							
Е	A S Y	A C C E S S							
	Abbeydale 760, 850, 950, 1050Edale MK2	Dovedale MK2							
S	TEP-IN	S H O W E R S							
	Berriedale	• Ferndale							



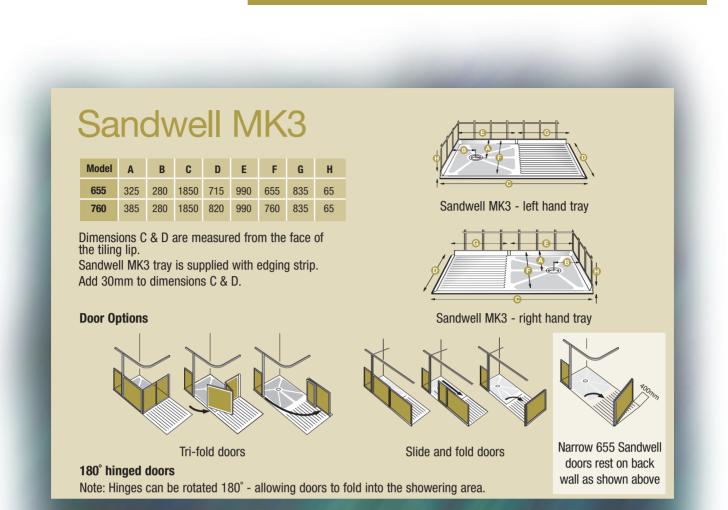

F

Westdale


Model	A	В	C	D	Е	F	G
655	325	285	1165	720	1100	655	65
760	385	285	1165	825	1100	760	65
1000	175	175	1065	1065	1000	1000	65

Dimensions C & D are measured from the face of the tiling lip.


Westdale tray is supplied with edging strip. Add 30mm to dimensions C & D.

Door Options

180° hinged doors

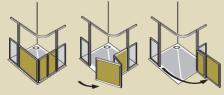


Tri-fold doors - Westdale 1000

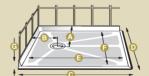

Westdale 655 - left hand tray



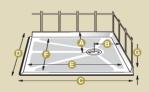

С

Е


S

Westdale 760 - right hand tray

Ν

0

Westdale 1000 - square

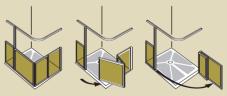





Note: Hinges can be rotated 180° - allowing doors to fold into the showering area on 655/760 models.

Height from tray to top of door 810mm Dimensions are in mm and for guidance only.



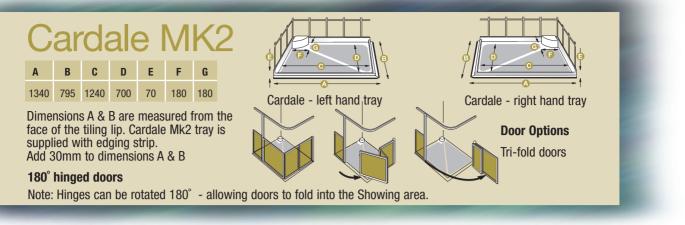

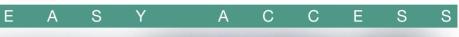


Abbeydale


Model	Α	В	C	D	Е	F	G
760	360	285	1200	860	1100	760	70
850	180	180	960	960	850	850	70
950	180	180	1080	1080	950	950	70
1050	180	180	1200	1200	1050	1050	70

Abbeydale tray is supplied with edging strip. Add 30mm to dimensions C & D. Dimensions C & D are measured from the face of the tiling lip


Door Options

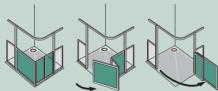


Tri-fold doors - double pole (850/950/1050 only)

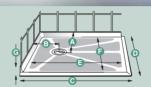


Abbeydale 760 - left hand tray

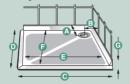

Abbeydale 950 - square tray

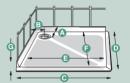


Abbeydale 850 - square tray

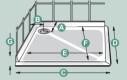

Abbeydale 1050 - square tray




Tri-fold doors - single pole (760 only)

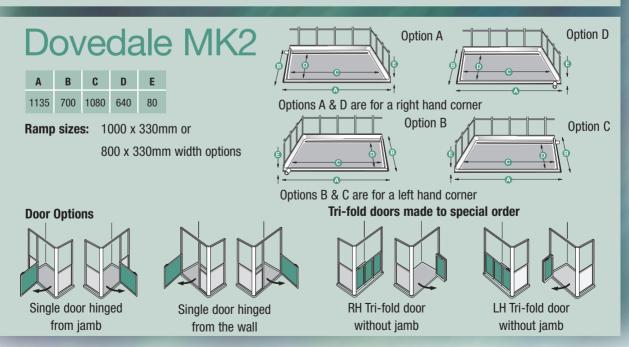

Height from tray to top of door 810mm Dimensions are in mm and for guidance only.



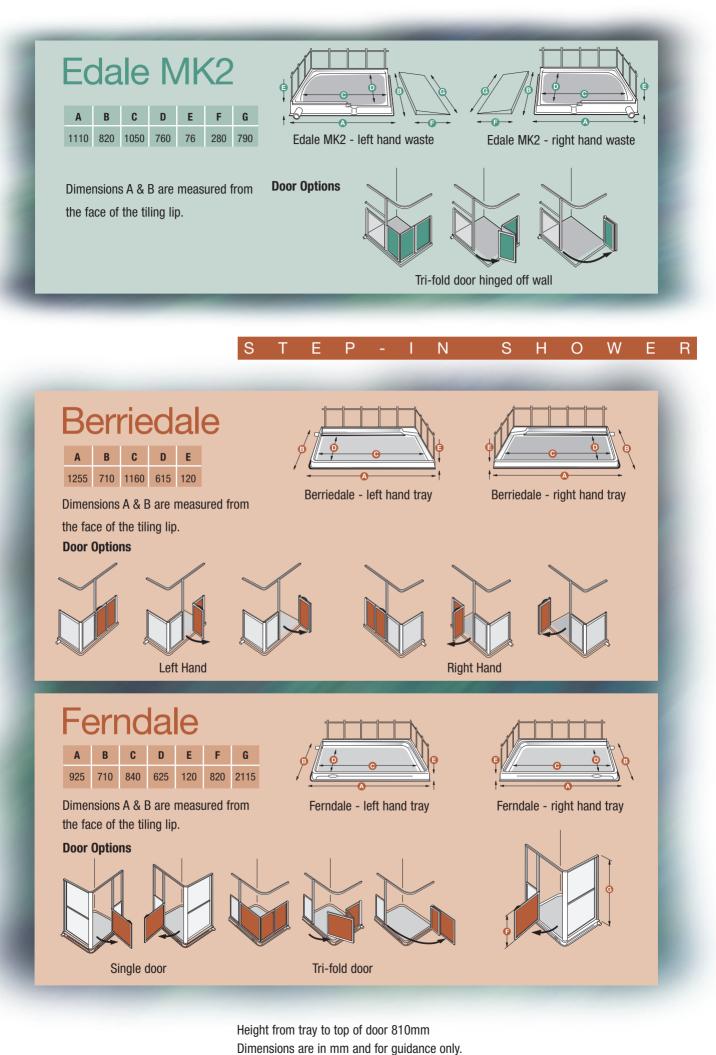

А

NFORMA

O N

Т

ECHNIC

Shower Chair Compatibility



Static Shower Chair WO2



Attendant Propelled Shower Chair W202AP/W204AP



Shower Chair

W201SP/S203SP

Polypropylene Shower Chair

W12MK2

To see the complete range in detail please ask for a copy of the Chiltern Invadex Shower Chairs brochure

Туре	Shower Tray	Shower Chair Compatibility								
Level Access	Sandwell 760	W02	W12MK2	W201SP	W202AP	W203SP (handrims removed)	W204AP			
Level Access	Sandwell 655	W02	W12MK2	W201SP (with rims removed)	W202AP		W204AP			
Level Access	Sandwell 760 with slide and fold doors	W02	W12MK2	W201SP	W202AP		W204AP			
Level Access	Sandwell 655 with slide and fold doors	W02	W12MK2	W201SP (with rims removed)	W202AP		W204AP			
Level Access	Westdale 655	W02	W12MK2		W202AP					
Level Access	Westdale 760	W02	W12MK2	W201SP	W202AP	W203SP (handrims removed)	W204AP			
Level Access	Westdale 1000	W02	W12MK2	W201SP	W202AP	W203SP (handrims removed)	W204AP			
Level Access	Cardale MK2	W02	W12MK2	W201SP (with rims removed)	W202AP		W204AP			
Easy Access	Abbeydale 760	W02	W12MK2	W201SP	W202AP	W203SP	W204AP			
Easy Access	Abbeydale 850	W02								
Easy Access	Abbeydale 950	W02	W12MK2	W201SP (with 90° footrests)	W202AP	W203SP (with 90° footrests)	W204AP			
Easy Access	Abbeydale 1050	W02	W12MK2	W201SP	W202AP	W203SP	W204AP			
Easy Access	Edale MK2	W02	W12MK2	W201SP (handrims removed)	W202AP		W204AP			
Easy Access	Dovedale MK2		NOT SUIT	ABLE FOR SHOWER	Chair Use					
Step-in	Berriedale	W02								
Step-in	Ferndale		NOT SUITA	BLE FOR SHOWER (CHAIR USE					





Solutions for living

Chiltern Invadex Ltd

126 Churchill Road Bicester, Oxfordshire OX26 4XD

Tel: 01869 365500 Fax: 01869 365588

E-mail: sales@chilterninvadex.co.uk Web: www.chilterninvadex.co.uk



Other products available from Chiltern Invadex include:

Showering and Bathing

Shower Cubicles Shower Accessories Shower Seats and Shower Chairs Shower Trolleys and Stretchers Walk-in and Easy Transfer Baths

Moving and Handling

Ceiling Mounted Hoists Wheelchairs Slings **Mobile Hoists Poolside Hoists**

Bicester is BS EN ISO 9001 approved (FM 40009). Bangor Northern Ireland is BS EN ISO 9002 approved (FM 21085).

Copyright © Chiltern Invadex Ltd

The design of the product shown on this brochure is vested in Chiltern Invadex Ltd and should not be copied or reproduced in anyway whatsoever without the express permission in writing of Chiltern Invadex Ltd.

As part of our commitment to new product development we reserve the right to alter specifications without prior notice. This is an original design of Chiltern Invadex Ltd in which Design Right subsists. Published by: Chiltern Invadex Ltd, 6 Wedgwood Road, Bicester, Oxon, OX26 4UL Registered in England No: 1182024 ISSUE 6 January 2003 In Northern Ireland: Chiltern Invadex (NI) Ltd, 5 Balloo Way, Balloo Industrial Estate, Bangor, County Down, BT19 7QZ Registered No: NI 19958