

CI/SfB	(72.36)	X	(U35)

Tip Up Shower Seat Range

Showering & Bathing Solutions



Made in the UK

TUW13BL shown with Sandwell tray and enclosure

Tip Up Shower Seat Range

Showering solutions for home, leisure and caring environments

The Chiltern Invadex range of tip up shower seats offer a safe and hygienic showering solution for those who are wheelchair dependent or have difficulty standing through showering.

The standard range comprises over 20 seat variations, including horseshoe and aperture designs to aid washing.

Manufactured with legs as standard, shower seats are available with optional armrests and backrest to provide additional comfort and security.

In addition to the standard range, bespoke seats can be produced to suit users with specific needs. This includes height, width and aperture alterations.

The robust design can safely withstand a maximum user weight of upto 31 stone (200kg). Subject to wall construction.

Standard Range



TUV18 Nylon reinforced PVC




TUV13BL Flat padded seat



Padded Horseshoe
TUV13LSH Small
TUV13LMH Medium
TUV13LLH Large



Padded Aperture
TUV13ALSAP Small
TUV13ALMAP Medium
TUV13ALLAP Large

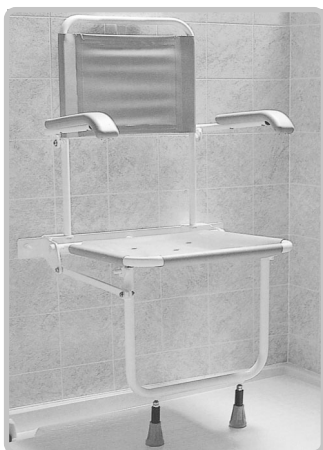

TUV19 Nylon reinforced PVC with backrest and folding arms



TUV14L Extra width padded bench 860mm
TUV15L Extra width padded bench 600mm
(horseshoe and aperture seat options available)



TUV19P Flat padded seat with backrest and folding arms (horseshoe and aperture seat options available)



The intuitive tip up design allows shower seats to fold neatly against the wall when not in use providing an ideal solution for multi-use bathing areas.



Shower seats are provided fully assembled and are easily fixed to virtually any wall construction.

An optional Surefix wall mounting kit (part number CHSFX) has been developed allowing safe and secure installation to timber studwork and soft blockwork.


It is important to remember that a wall mounted shower seat is only as secure as the wall mountings and the wall to which it will be fixed. We recommend that if there is any doubt a reputable builder/installer is consulted.

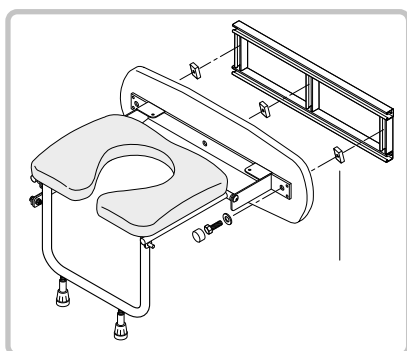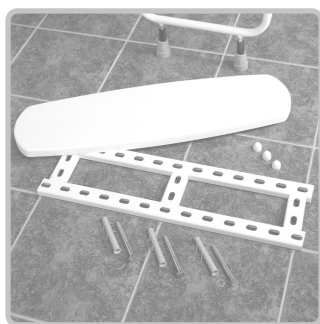

Each seat is supplied with non-slip adjustable feet which can be used to vary the seat height (at the time of fixing).

Standard models are supplied with legs which can be fitted to provide a seat height of between 458mm and 509mm.

At the time of order legs can be requested which allow height adjustment of between either 407mm to 458mm or 509mm to 560mm.

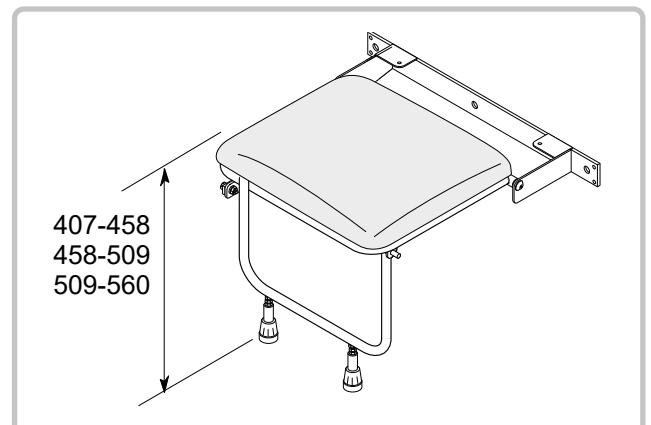

Adjustable feet can also be used to stabilise the seat on uneven flooring.

Safety

Designed and tested with user safety in mind, the Chiltern Invadex range of shower seats meets the robust safety requirements of BS EN 12182È

TUW18 shown with wet floor installation




Solutions for Living

Chiltern Invadex (UK) Ltd
126 Churchill Road
Bicester
Oxfordshire
OX26 4XD
Tel: 01869 365500
Fax: 01869 365588

E-mail: sales@chilterninvadex.co.uk
Web: www.chilterninvadex.co.uk

Overseas enquiries
Tel: 0044 1869 365500
Fax: 0044 1869 365588

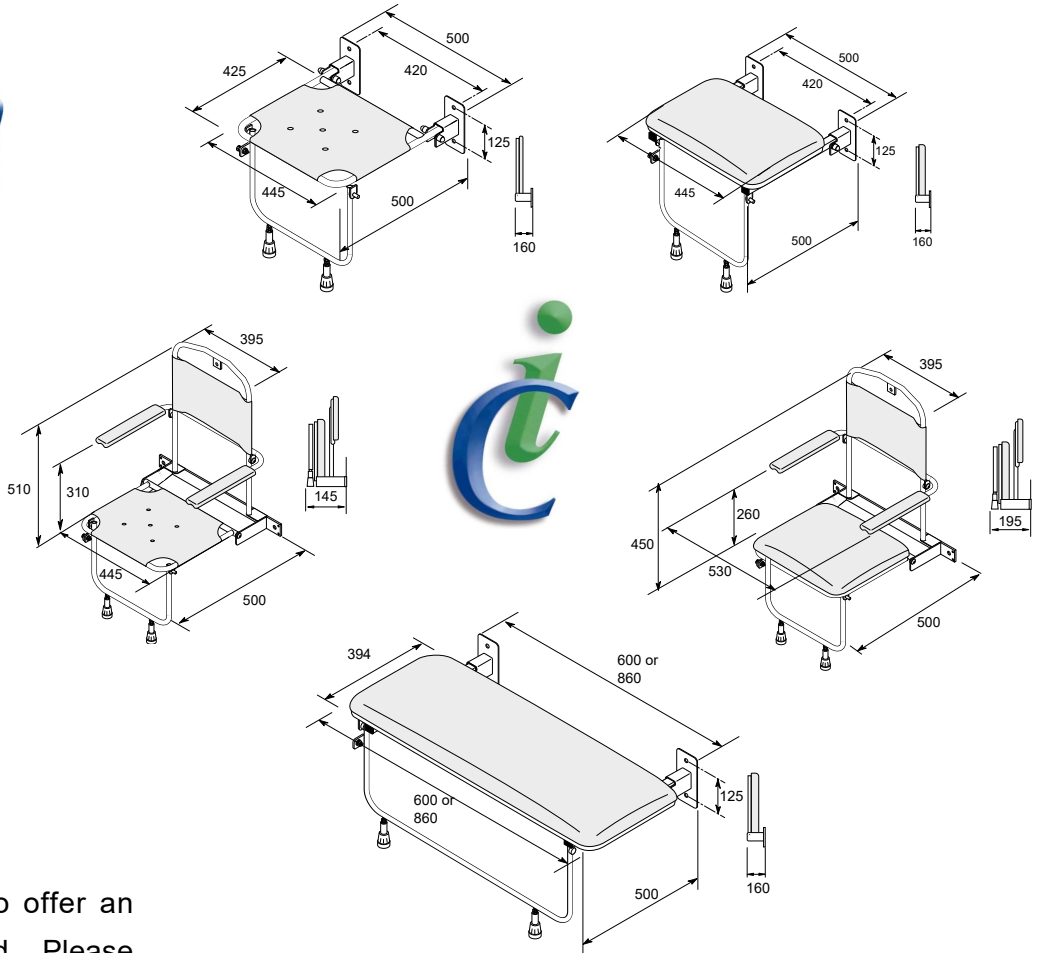


Installation

Chiltern Invadex are able to offer an installation service if required. Please contact our service department for a quotation on 01869 365500 Opt 3.



A range of shower enclosures, cubicles and shower commode chairs are also available, to obtain a brochure please contact customer services on 01869 365500 Opt 1 or email sales@chilterninvadex.co.uk.

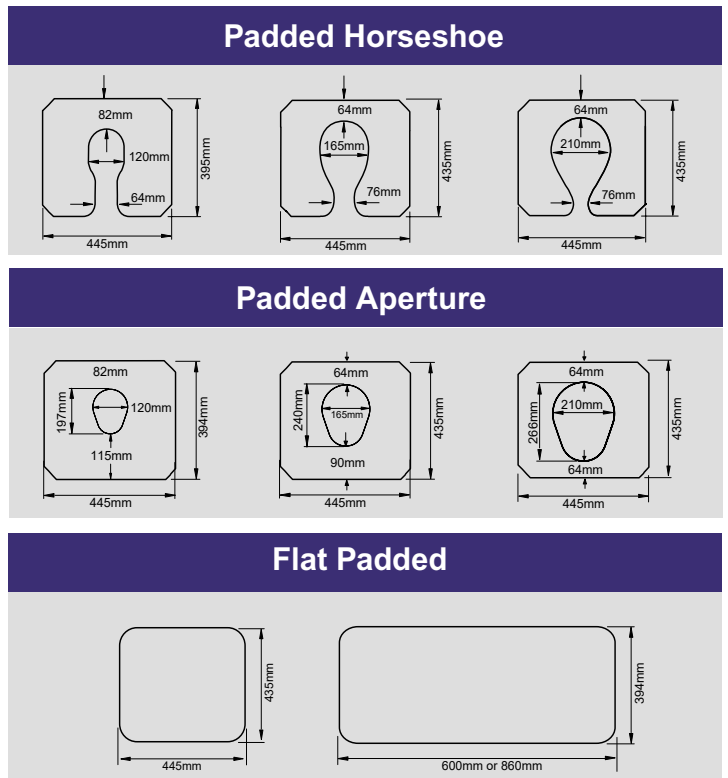

FM 553189
BS EN ISO 9001:2008